



TIBIA INTRAOSSEOUS (EZ-IO®)


Continued...

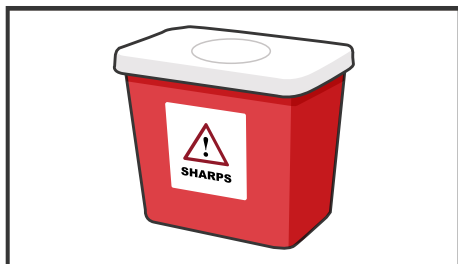


- 13** **PLACE** the stylet in a sharps container.

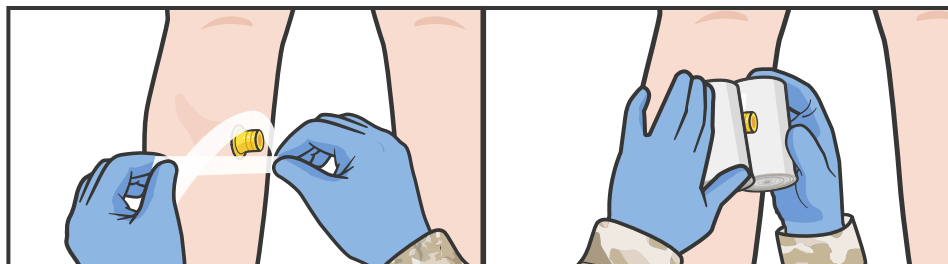


- 14** **SECURE** the site with an EZ-IO stabilizer, if available.
NOTE: If a stabilizer is unavailable, secure with bulky dressing and tape to prevent elevated extension tubing from becoming dislodged during casualty care or movement.

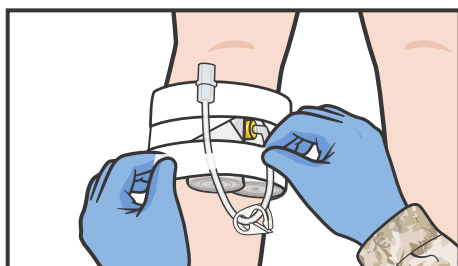


- 15** **CONNECT** primed extension set through the 90-degree IO tubing to the Luer lock EZ-IO hub.

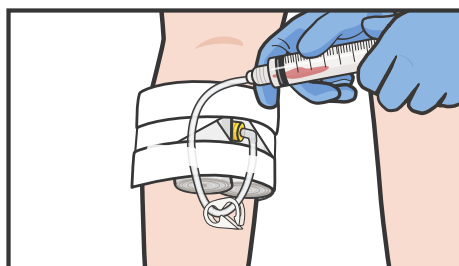


- 16** **ASPIRATE** to confirm needle is in the marrow (should see flash of blood-tinged material) and flush catheter with 10 ml flush using a rapid infusion.

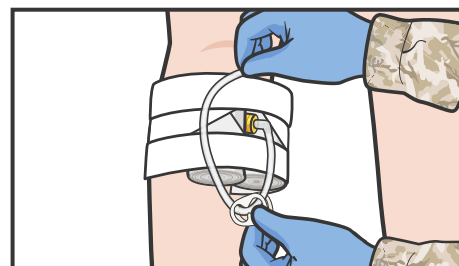

- 17** **ASSESS** for signs of infiltration or complications.

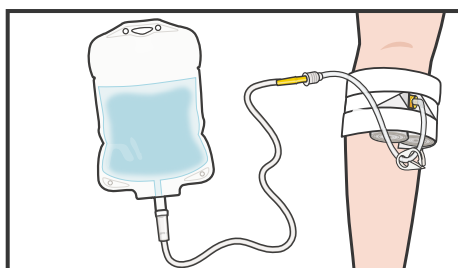

- 18** **PREPARE** for fluid or medication administration.

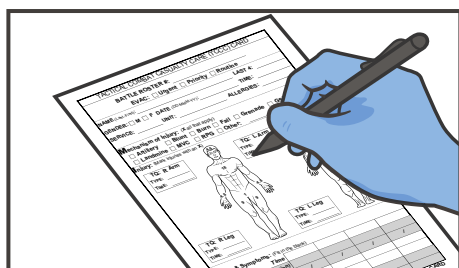

- 19** **DOCUMENT** all findings and treatments on a DD Form 1380 TCCC Casualty Card and attach it to the casualty.