



STERNAL INTRAOSSEOUS (FAST1® IO)

Continued...

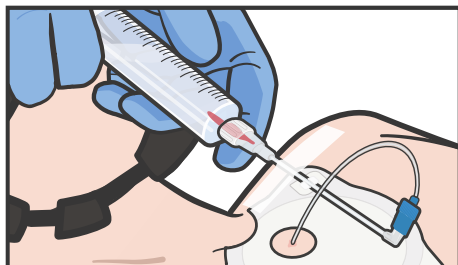

13 **CONNECT** the syringe, aspirate to confirm needle is in the marrow (should see flash of blood-tinged material) and flush the infusion tube with 5 ml of sterile IV solution.


NOTE: This is necessary to flush the bone plug or any obstructing tissue.

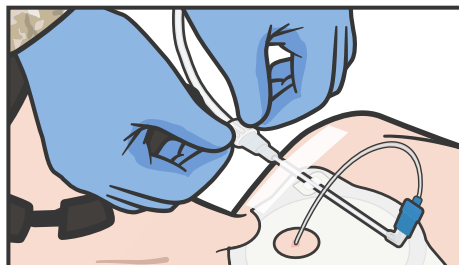


14 **CONNECT** the IV infusion tubing and initiate flow of fluids.



15 **ATTACH** the protective dome to the target patch and secure with tape.

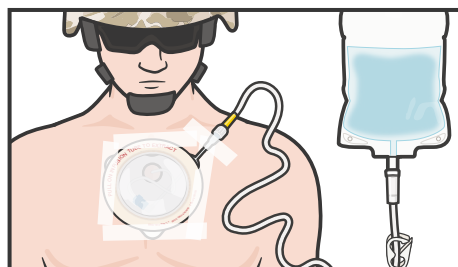

16 **ADMINISTER** IV fluids and/or medications as needed.

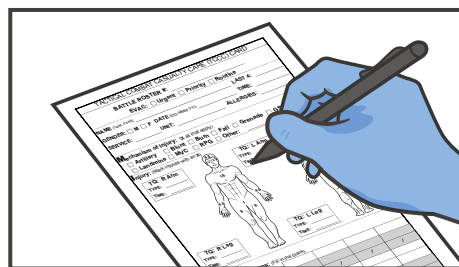

17 **DOCUMENT** all findings and treatments on a DD Form 1380 TCCC Casualty Card and attach it to the casualty.